

Sonicaid™
Fetal Dopplers and Monitors

HUNTLEIGH



...with people in mind

sonicaid™ Fetal Monitoring


Performance... for life

Performance by Sonicaid


Commitment to Excellence and Durability

Performance for life is at the core of every Sonicaid product. While our products are primarily designed and manufactured to assist in the safe delivery of a new life, absolute care and attention to detail is paramount to ensure that the high levels of performance are consistent throughout the life of the product.

Manufactured in the
United Kingdom...

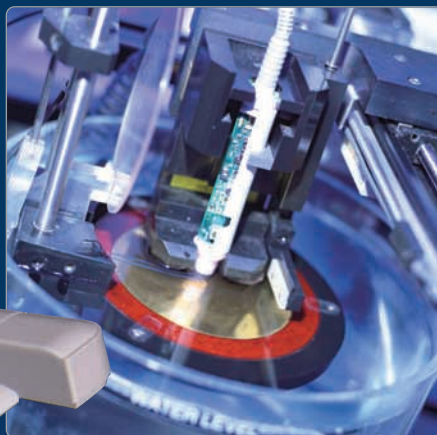

Designed in the
United Kingdom...

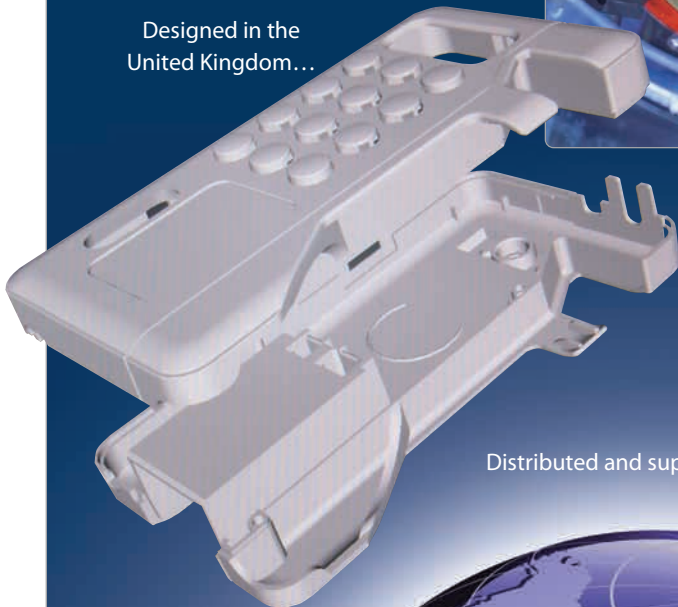

Distributed and supported globally.



sonicaid Dopplers

HANDHELD 2MHz & 3MHz DOPPLERS



The latest generation of the world renowned hand held fetal Dopplers offer even greater performance, quality and value for money.

The ultimate in obstetric pocket Doppler for obstetricians, midwives, doctors and other specialists.

The original...
often copied...
never equalled.



Approved accessories and consumables available – ask about our contract prices

waterproof

Range includes waterproof models



sonicaid **BD4000xs**

ENTRY LEVEL FETAL MONITORING



The Sonicaid BD4000xs is designed to deliver cost-effective monitoring without compromising quality, precision, performance and ergonomics. Ideally suited for the antenatal monitoring environment.

Twins capable as standard on wide format paper.



Approved accessories and consumables available – ask about our contract prices



sonicaid **TEAM**

ADVANCED LEVEL MONITORING WITH FHR ANALYSIS



The Sonicaid Team series combines standard CTG monitoring with the world renowned Dawes/Redman antepartum CTG analysis, providing an objective, numeric and consistent assessment of the CTG.

Ideally suited for antenatal screening, giving clinicians advanced warning of fetal compromise.

Twins, tele-fetal monitoring and intrapartum variants available.



Approved accessories and consumables available – ask about our contract prices



CTG Analysis



Tele-Fetal Monitoring option available

sonicaid **FM800** *Encore*

ACUTE FETAL/MATERNAL INTRAPARTUM MONITORING



The Sonicaid FM800 Encore series is designed to provide flexible solutions for high-risk fetal and maternal assessment.

With integrated maternal vital signs monitoring (FM830E), the FM800 Encore series offers enhanced care throughout labour and delivery.

Includes antepartum CTG analysis and Intrapartum Trend functionality.



Approved accessories and consumables available – ask about our contract prices



CTG Analysis



Trend report function for enhanced labour management



Wireless transducer option available

sonicaid Dopplers

Often copied... never equalled

SONICAID FETAL MONITORING

Sonicaid Dopplers - the ultimate in Obstetric pocket Dopplers for obstetricians, midwives, doctors and other specialists.

Sonicaid Value

The exceptional build standard and quality of a Sonicaid Doppler not only gives outstanding performance, but also durability and a longer product life.

Our high sensitivity probes ensure optimum clinical performance even with difficult assessments such as early gestation or raised BMIs.

Combine this with the added advantage that a Sonicaid Doppler will often outlive many of its competitors – You may be pleasantly surprised to find that the professional's choice of fetal Doppler offers unrivalled value.

*Don't just buy a Doppler
...invest in a Sonicaid*



Sonicaid Fetal Dopplers

The best value fetal Doppler?*

- Sonicaid launched the world's first handheld fetal heart detector in 1971
- Huntleigh launched the world's first rate display Doppler in 1990
- Today, Huntleigh Sonicaid leads the world in Doppler innovation, sensitivity, quality and reliability

* Based on total life cost.

sonicaid One

RATE DISPLAY DOPPLER WITH FIXED 2MHz PROBE



This compact fetal Doppler has been designed to fit easily into the clinician's pocket. Its wide beam coverage enables quick location of the fetal heart from early gestation to full term.

- Innovative and stylish
- Integrated probe and cable storage
- Low running costs, excellent battery life
- Pocket ergonomics
- Large fetal heart rate display

sonicaid D920 and D930

AUDIO DOPPLERS WITH FIXED 2MHz OR 3MHz WATERPROOF PROBE



The D920, with its high sensitivity probe and excellent audio quality, ensures easy and reliable detection throughout pregnancy. Our latest D930 model meets the growing demand for a high sensitivity 3MHz Doppler for sensitive detection in early gestation.

- High sensitivity probes for improved ease of use
- Waterproof probes suitable for waterbirth
- Easy clean probe for improved infection control



The Sonicaid HS range offers an unmatched combination of unique automatic noise reduction, acclaimed fetal heart sounds and unequalled sensitivity. Allows users to more easily detect fetal sounds in early gestation and mums with raised BMIs (Body Mass Index) than ever before.

sonicaid FDI and FD3

RATE DISPLAY DOPPLERS WITH FIXED 2MHz OR 3MHz WATERPROOF PROBE MODELS



The FD1 combines all the powerful features of the D920 with the added advantage of displaying fetal heart rate. The high sensitivity 2MHz probe provides reliable detection throughout pregnancy. Our latest FD3 model meets the growing demand for a high sensitivity 3MHz rate display Doppler for sensitive detection in early gestation.

- Advanced fetal heart rate (FHR) processing
- Active noise reduction for crystal clear fetal heart sounds
- Waterproof probes suitable for waterbirth use
- Advanced power management to conserve battery life
- Easy clean probe for improved infection control
- Long robust coiled cable enables optimum transducer position in all circumstances



waterproof
probe as standard

sonicaid FD2

RATE DISPLAY DOPPLER WITH INTERCHANGEABLE OBSTETRIC AND VASCULAR PROBE OPTIONS



This top of the range handheld rate display Doppler features high sensitivity interchangeable probes. 2MHz for general obstetric use, 3MHz ideally suited for detecting the fetal heart in early gestation. A unique interface to our portable Printa II provides an FHR trace print-out, complete with fetal movement marker.

- Probe connector for interchangeable probe options, including the vascular range
- 3 FHR processing modes for optimum performance in the widest range of signal conditions
- Real time fetal heart rate data output
- Easy clean probe for improved infection control



sonicaid MD200

RATE DISPLAY DOPPLER WITH INTERCHANGEABLE OBSTETRIC AND VASCULAR PROBE OPTIONS



The MD200 is an innovative solution for clinicians requiring a robust versatile and rechargeable desk top Doppler unit. Functionally equivalent to the Sonicaid FD2, features also include:

- Mains/Rechargeable battery operation
- Two loudspeakers providing superior audio
- Probe and cable storage compartment



Integral probe parking and protective bumper protects the high sensitivity probe from damage.



High quality speaker gives high quality audio.

sonicaid™ BD4000xs

...simplicity

ENTRY LEVEL FETAL MONITORING

The BD4000xs fetal monitor is ideal for antenatal clinic and doctor's office use.

Its compact size and styling belie its ruggedness and robust design, using the same materials, construction and design expertise as our world-renowned handheld Dopplers.

- Twins capability as standard
- Unique no compromise full size twin FHR scale printing
- User adjustable, tachycardia, bradycardia, loss of contact and cross – channel rate verification alerts
- Software System Interface
- User-customisable trace annotation facility
- Automatic fetal movement detection
- Intuitive user interface and controls simplifying training
- Maternally sensed fetal movement marker
- Extra long plain thermal paper packs reduces running costs and guarantees 100% trace registration accuracy.
- Trolley and wall mount options



Approved accessories and consumables available – ask about our contract prices



Connects to Sonicaid Centrale and other CTG viewing and archiving software systems. Interface included as standard at no extra cost.

NO COMPROMISE TRACES

"Dual channel" A4 width Z-fold chart printout for unique optimum twins presentation on plain thermal paper for low cost, high accuracy printing (standard width format paper can also be used).

INTUITIVE USER INTERACTIVE TEXT DISPLAY

A feature packed fetal monitor, including trace annotation facility, user configurable alarms, automatic fetal movement detection, and external communication ports all included as standard.

ERGONOMIC DESIGN

The BD4000xs CTG is the monitor of choice, whilst being both robust and a remarkable performer, it is also an elegant and slim looking system.

SINGLE OR TWINS CAPABILITY

Clear visual displays of FHR, contractions data and signal quality indicators.

EASY TO USE

Only one button press to start monitoring.



sonicaid™ TEAM

...CTG analysis

ADVANCED LEVEL FETAL MONITORING WITH CTG ANALYSIS

The Sonicaid Team range is designed for advanced fetal monitoring. Developed in collaboration with leading obstetric specialists, the small, easy to use and lightweight unit can be transported effortlessly; within hospitals, to clinics, doctor's offices or even to the home.

- Twin fetal heart rate scales provide clear differentiation between two heart rate traces
- Software systems interface as standard
- Six hours of trace storage
- Maternally sensed and automatic (Actogram) fetal movement detection
- Alarms
- Twins (TeamCare Duo), Tele-Monitoring and intrapartum (Team IP Trend) variants
- TREND provides clear data on periodic changes over time in the FHR. For use during the 1st stage of labour. TREND highlights changes in baseline rate, Short Term Variation (STV) and decelerations



Approved accessories and consumables available – ask about our contract prices



CTG Analysis



Tele-Fetal Monitoring option available



Connects to Sonicaid Centrale and other CTG viewing and archiving software systems. Interface included as standard at no extra cost.

Change recorder paper easily and quickly.

Flat bed recorder design with large 'write on' surface.

Large illuminated display provides clear monitoring at a glance. Scrolling display of fetal heart rate and uterine activity traces.

Records can be annotated from on-board library of clinical notes.


Easy to use menu system allows fast access to important function.

Toco zero at the touch of a button.

Small, compact design for easy positioning.

New sealed high sensitivity toco transducer.

Colour coded transducers for easy recognition.

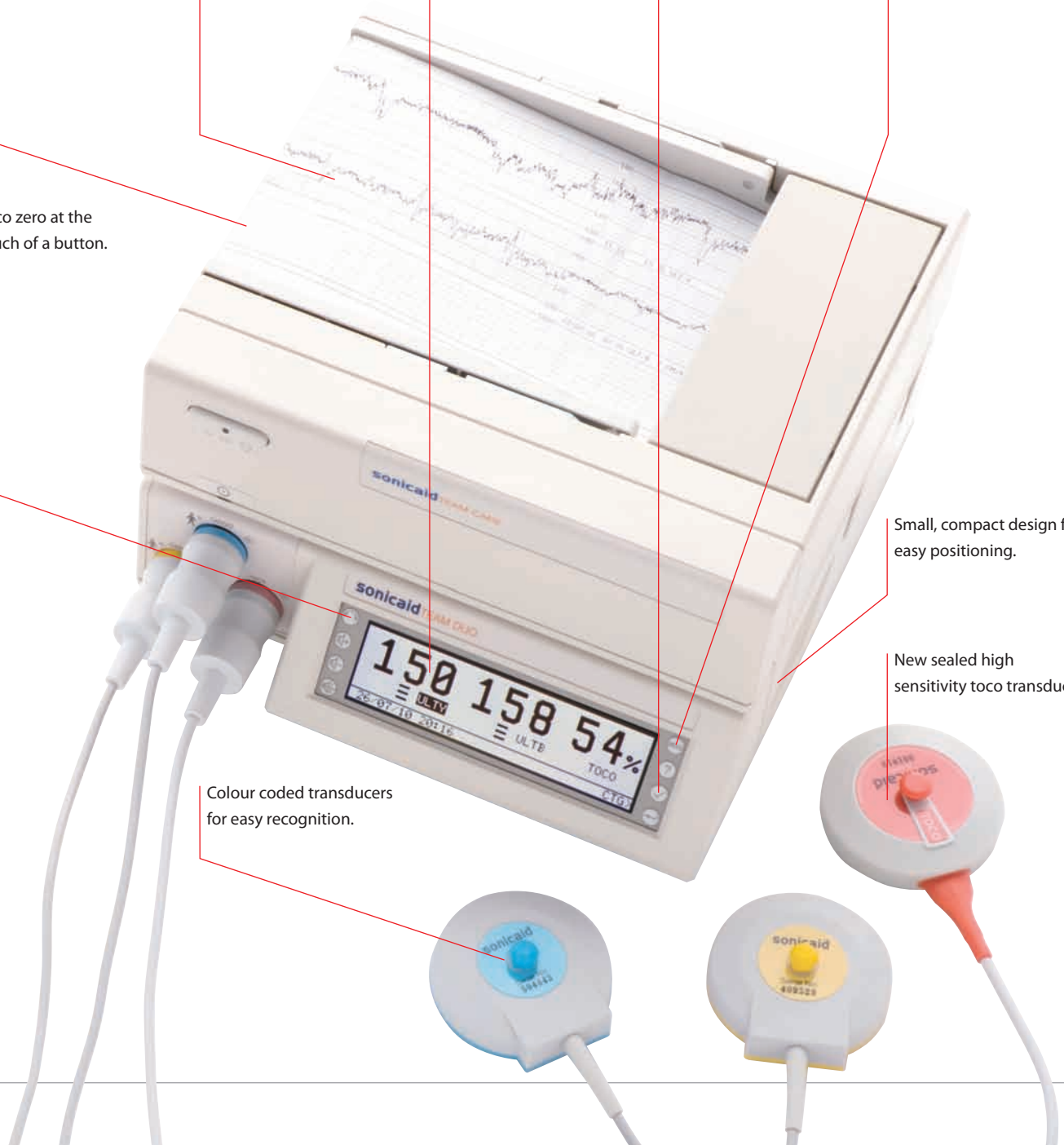

sonicaid™ FM800 Encore

...great performance

HIGH PERFORMANCE ACUTE FETAL/MATERNAL INTRAPARTUM MONITORING

The Sonicaid FM800 Encore series high performance fetal and maternal monitors, uses today's advanced technology to produce the most accurate acquisition of data. The Sonicaid FM800 Encore series gives you accuracy and clarity in every situation. In the fast-paced environment of the labour ward, you need equipment that's easy to use and data that is quick to read.

- Using the very latest in ultrasound technology, the Sonicaid FM800 Encore gives increased performance. Combined with our 'Locate & Track' technology, the FM800E makes it easy to locate and maintain contact with the fetal heart in mothers of all sizes
- Superbly accurate fetal ECG that works in the most demanding of environments – perfect for the busy labour room
- User defined alarm settings
- Unique and powerful Antepartum CTG Care analysis provides an objective, numerical and consistent assessment of the CTG
- TREND intrapartum function provides clear data on periodic changes over time in the FHR. For use during the first stage of labour
- Optional wireless transducer system available (Sonicaid Freedom)



Approved accessories and consumables available – ask about our contract prices



CTG Analysis



Trend report function for enhanced labour management



Wireless transducer option available



Connects to Sonicaid Centrale and other CTG viewing and archiving software systems. Interface included as standard at no extra cost.

Twin ultrasound plus separate fetal/maternal ECG channels come as standard. Both fetal and maternal ECG waveforms can be displayed on screen and as a "snapshot" printout for improved hook-up and monitoring.


Connecting to Sonicaid Centrale (now also via LAN) enables electronic data storage, reducing paper print-outs, giving even lower running costs and less environmental impact.

Clear colour screen graphics make it easy to monitor both mother and baby's vital signs at a glance. Transducers match on-screen colours for fast, error-free recognition.

Ergonomic design creates a soft and unobtrusive presence, while the slim-line trolley is easy to manoeuvre, with integrated cable management and paper catching tray.

Easy to learn and use. One-touch operation gives instant access to alarms and critical settings.

Quick-view trace display allows instant review at the touch of a button.

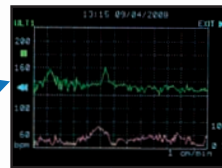

Records can be annotated from Easinotes, a programmable library of clinical notes, available at the touch of a button.

User selectable twins trace format: -20 beat trace offset or separate FHR scales

Memory feature enables paper change with no loss of data. Compatible with other brands of paper.

Toco zero at the touch of a button combined with our unique auto-zero function.

Colour coded transducers coordinate with screen for fast, error-free recognition. Waterproof design makes for thorough and easy cleaning – essential for infection control.

sonicaid™ Obstetric

Accessories and consumables ... maximising fetal monitoring performance

Huntleigh's medical products are renowned for reliability and superior quality, values which are reflected in our accessory and consumables product range, backed by experience and support from one of the world's leading healthcare companies. Whatever your needs, Huntleigh will work with you to find the optimum accessories and consumables that are right for your applications, helping you seamlessly manage important fetal monitoring requirements. For the very best value, consumable supply contracts are available in selected markets.

JUST A FEW EXAMPLES OF OUR ACCESSORIES INCLUDE:



Gel – it is essential to use the highest quality ultrasound gel to get the best performance out of your Doppler or fetal monitor. Our gel is available in a range of packs, from 60ml tubes to 5ltr refill containers.



Our "Longlife" CTG paper is guaranteed to retain clarity of CTGs for up to 25 years. Our 45m long packs last 3 x times longer than other makes, giving 72 hours continuous trace recording (at 1cm/min). For use in Team and FM800/FM800E series fetal monitors only.



Huntleigh's transducer belts prove their value every day both in labour and delivery rooms. Options include latex free reusable or disposable belts with stretch button holes for quick and easy transducer attachment.

Huntleigh's comprehensive accessories and consumables catalogue is available to download at:

www.huntleigh-diagnostics.com

or by request via:

sales@huntleigh-diagnostics.co.uk

By choosing our range of accessories and consumables, you ensure that you will get the very best performance and longevity from your product.

Note that use of unauthorised 3rd party accessories may effect the performance, quality and reliability of your product and may invalidate your warranty.

HUNTLEIGH ...with people in mind

Huntleigh Healthcare Ltd. - Diagnostic Products Division

35 Portmanmoor Road, Cardiff, CF24 5HN, United Kingdom

T: +44 (0)29 20485885 F: +44 (0)29 20492520 E: sales@huntleigh-diagnostics.co.uk

W: www.huntleigh-diagnostics.com



Registered No: 942245 England & Wales. Registered Office:
ArjoHuntleigh House, Houghton Hall Business Park, Houghton Regis, Bedfordshire, LU5 5XF
©Huntleigh Healthcare Limited 2012

AN ARJOHUNTLEIGH COMPANY, MEMBER OF GETINGE GROUP

© and ™ are trademarks of Huntleigh Technology Limited
As our policy is one of continuous improvement, we reserve the right to modify designs without prior notice.

LIT Number: 709419-4